


Wound Care Line

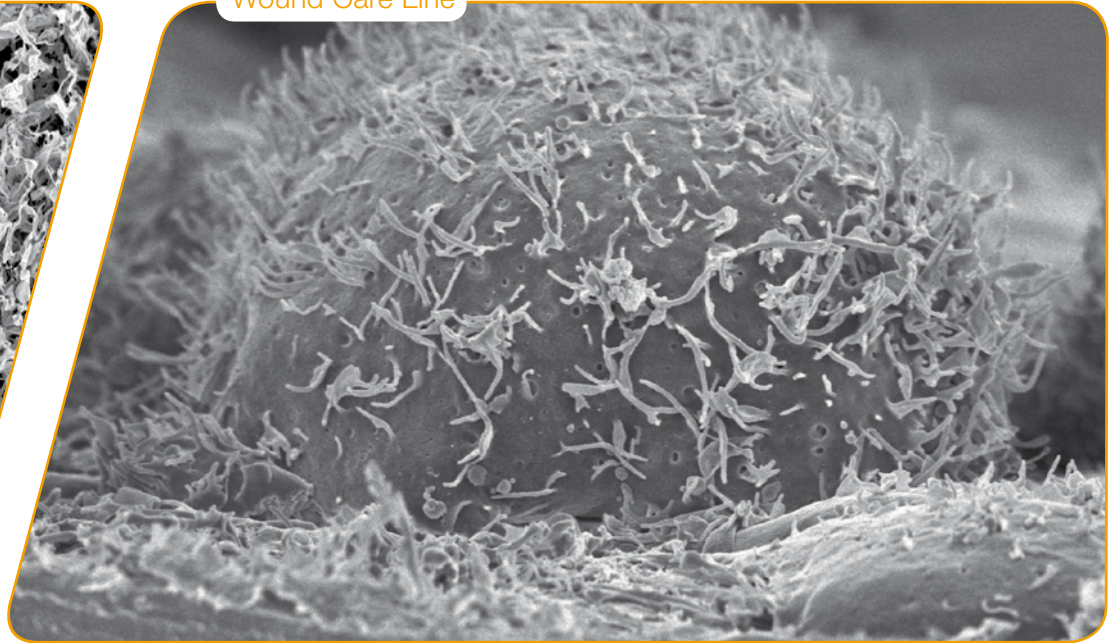

MB-Collagen

Collagen dressing for wound care

Properties:

- Highly effective local hemostatic agent.
- Promotes wound healing.
- Fully resorbable within 1-2 weeks.
- Reduces the number of dressing changes necessary.
- Easy and simple handling.

MB-Collagen

CE 0499



MB-Collagen is a collagen based dressing for wound care of traumatic and slowly healing wounds like Decubitus or Ulcus cruris.

The effectiveness of MB-Collagen is based on one hand on the chemotactic effects of the collagen matrix leading to a settlement of fibroblasts on the wound bed. On the other hand exogen collagen promotes binding of proteases. Both mechanisms are acting complementarily thus providing the basis for efficient wound care. MB-Collagen is also suitable for combination with other components like e.g. bactericidal agents.

Further Properties:

- Class III medical device according to European guideline 93/42/EWG.
- Certified according to DIN EN ISO 13485:2012.
- Production following good manufacturing practice guidelines (GMP).
- Outstanding biocompatibility.
- Gamma sterilized.

Specifications

Art-Nr.	Size	Contents
REF 0505C	50x50x6 mm	3 pieces per box
REF 0510C	50x100x6 mm	3 pieces per box
REF 1010C	100x100x6 mm	3 pieces per box

